

GRACIE

LINET

The Next Level of Patient
and Gynecologist Experience



Modern Gynecological Workplace

GRACIE

A premium gynecological solution designed to ease patient tension and increase clinical efficiency in workplace.

Optimal Workplace

Easy switching between individual settings for up to three doctors. The doctor's office can simply be shared.



Ultrasound examinations can be performed on the chair without the need to move the client to another location.



The instrument trolley with ergonomic storage spaces is easy to access and move.



Modern Equipment



Integrated colposcope securely stored under the seat and ready for immediate use.



Smooth Operation



Position Selection with a wireless controller and subsequent adjustments with foot and hand controls.



Quick positioning including pre-set buttons for the most frequent chair adjustments.



Easy access to the chair even for less mobile woman thanks to complete removal of the footrests.



FULL HD LCD and wireless transfer of photos and videos to the doctor's PC.

TECHNICAL PARAMETERS

Height adjustment


Access height (seat)	520 ±5 mm
Max. height (seat)	1,150 ±5 mm

Maximum load

Patient load	180 kg
Maximum load bearing capacity of the chair	210 kg

Dimensions

Chair width – seat	630 ±5 mm
Chair width – total	760 ±5 mm
Chair length – max.	1,780 ±5 mm

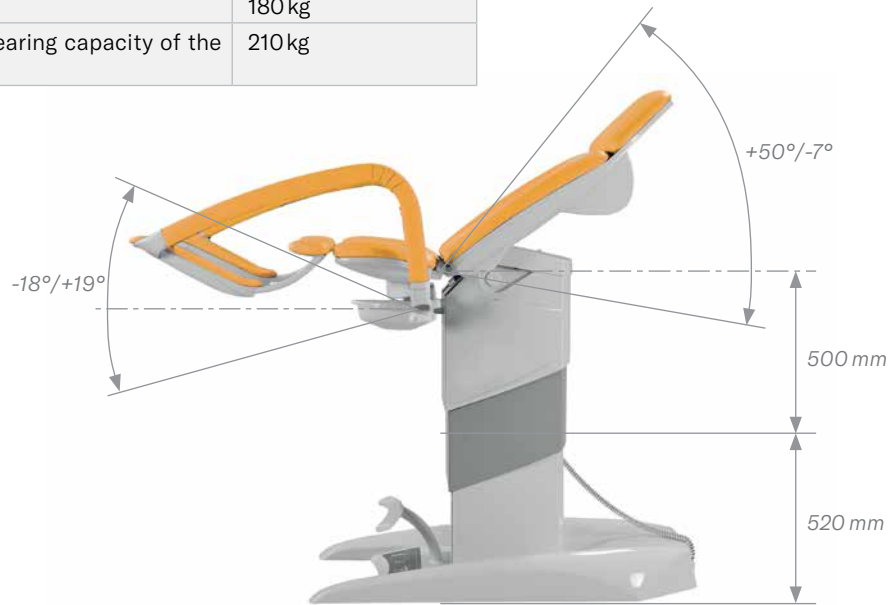

Electrical parameters

Voltage	100 –230 V
Frequency	50/60 Hz
Applied parts type	CF
Protection class	I
Degree of protection	IP X4

Power adjustment

Back section adjustment angle	-7°/+50°
Foot rest adjustment angle	-18°/+19°
Seat section adjustment angle	+5°/+10°



Members of LINET Group

LINET spol. s r.o.

Želevčice 5 | 274 01 Slaný | Czech Republic

tel.: +420 312 576 400 | fax: +420 312 522 668 | e-mail: info@linet.com | www.linet.com



www.linet.com